

Medimax System

For Automated,
Tissue Dissociation

STANDARDIZED
AND EFFICIENT
SAMPLE PREPARATION
TO OBTAIN VIABLE
CELL SUSPENSION
FROM SOLID TISSUE



code
79400



code
79300

code
79500



www.syntec-international.com

Medimax System

For Automated, Tissue Dissociation

Standardization and Increased Efficiency in Sample Preparation

A small piece of tissue (up to 10 mm³), free of fat and necrosis, is inserted into the Medicons together with **the correct amount of liquid** (~1.0 mL for codes 79300 and 79400; ~16.0 mL for code 79500) of the desired suspension buffer (for example, PBS)

The Medicons is then inserted into the Medimax and the machine is started. The machine can be run from 10 seconds to 4 minutes depending upon the tissue type and the cell suspension you want. As the tissue is disaggregated, the cells pass through the screen into the suspension liquid in the bottom of the Medicons. Once the tissue is dissociated, the Medicons is removed from the Medimax. A needleless syringe is then inserted in the syringe port to recover the cell suspension. The Medicons may be rinsed with suspension buffer to ensure maximum cell recovery. Select the appropriate Filcons (or Filcons-Max if a 50 mL tube is used) pore size for isolation of the desired cells.






Medicons

Medicons are disposable polyethylene chambers containing an immobile stainless steel screen with approximately 100 hexagonal holes. Around each hole are six microblades designed for efficient cutting of hard and soft tissues.



Filcons

After disaggregation of the tissue in the Medicons, the cell suspension is filtered using a Filcon disposable filter device. Depending upon the chosen Filcon pore size, cells or cell aggregates pass through, while larger particles are trapped. There are two types of Filcons: a syringe-fitting form which can be used with either suction or compression, which works well for quantities over 5 mL, and a cup type that works well for small quantities.

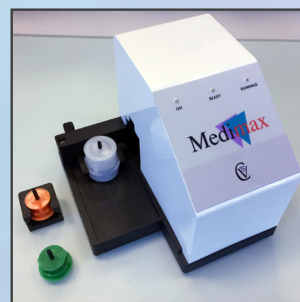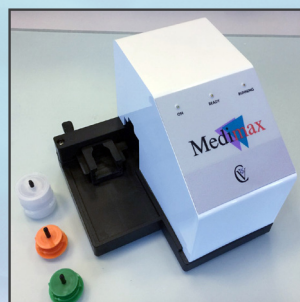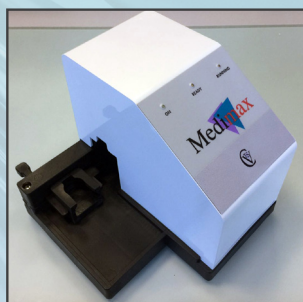

TECHNICAL TIPS:

MEDICONS

- Never place the tissue into a dry Medicons.
- Do not overfill the Medicons with suspension buffer. This will cause the tissue to float and can result in inadequate disaggregation: dilution must be between 1 and 1,2 mL for standard Medicons (codes 79300 and 79400) and between 15 and 16 mL for Large Medicons (code 79500).
- If the amount of tissue is very small, the suspension buffer can be replaced with staining solution for maximum cell recovery.
- As the rotating element passes over the microblades, it tends to blunt the edges, reducing the cutting efficiency.
- Do not reuse the Medicons for successive samples.
- Large Medicons (big white) is available in 50 µm mesh size, while Medicons standard (green and orange) are available (sterile or not sterile) in two mesh sizes:

35 µm	Recovery of cell sizes up to 160 µm
50 µm	Recovery of cell sizes up to 230 µm

MEDIMAX

The Medimax system is a sample preparation system for the automated, mechanical disaggregation of solid animal or plant tissues for flow cytometric analysis, cell culture, or DNA amplification techniques.

The Medimax system makes tissue preparation safe and independent of operator technique, and facilitates standardizing the preparation of all tissue types. Times and protocols can even be standardized across labs.

The Medimax system consists of three components:

- Medimax
- Medicons
- Filcons & Filcons Max

The compact Medimax unit is the heart of the system. It works with all types of Medicons and operates at a constant speed of approximately 80 rpm.



REG. N. 623
UNI EN ISO 9001:2008

Sede Operativa Via Marco Polo 28 - 10090 Bruino (TO) Italy
info@ctsv.biz - Tel. +39 011 9049381 - Fax +39 011 9055288 - www.ctsv.biz

